



## Product Data Sheet

### **DuPont™ Duolite™ AP143/1073 Resin**

Pharmaceutical Grade Anion Exchange Resin Powder (Cholestyramine Resin USP)

#### **Description**

DuPont™ Duolite™ AP143/1073[1] resin is an insoluble, strongly basic, anion exchange resin in the chloride form supplied as a dry, fine powder. Duolite™

AP143/1073 Resin is suitable for pharmaceutical applications either as an active ingredient or as a carrier for acidic (anionic) drug substances. A monograph for Cholestyramine Resin USP appears in the United States Pharmacopoeia/National Formulary.

Duolite™ AP143/1073 resin conforms to the compendial specifications.

A Drug Master file for this product is maintained with the United States Food and Drug Administration.

Duolite™ AP143/1073 Resin can be identified by infrared spectroscopy, as shown in the example Figure 1. Letters of authorization granting access to the file by FDA in support of NDA and ANDA submissions will be provided upon request.

Duolite™ AP143/1073 Resin is manufactured in accordance with Good Manufacturing Practices (cGMP) for bulk pharmaceutical chemicals.

<sup>1</sup> The use of DuPont™ AmberLite™ and DuPont™ Duolite™ pharmaceutical grade ion exchange resins as components of drug formulations is subject to the Food, Drug, and Cosmetic Act as amended.

## Typical Properties

Physical Properties	
Copolymer	Styrene-divinylbenzene
Type	Strong base anion
Functional Group	Quaternary amine
Physical Form	White to buff-colored, fine powder
Chemical Properties	
Ionic Form as Shipped	Cl <sup>-</sup>
Loss on drying <sup>1</sup>	≤ 12.0%
Identity (by IR spectrum) <sup>1</sup>	Identical to USP reference standard
Heavy metals <sup>1,2</sup>	≤ 0.002%
pH of slurry <sup>1</sup>	4.1–6.0
Residue on ignition <sup>1</sup>	≤ 0.1%
Dialyzable quaternary amine <sup>1</sup>	0.05%
Chloride content <sup>1</sup>	13.0–17.0%
Sodium glycocholate exchange capacity <sup>1</sup>	1.8–2.2 g/g
Trimethylamine <sup>1</sup>	≤ 20 ppm maximum
Organic volatile impurities <467> <sup>1</sup>	Meets requirements
Particle Size <sup>§</sup>	
< 150 μm <sup>1</sup>	≥ 95%
< 75 μm <sup>1</sup>	≥ 65%
< 45 μm <sup>1</sup>	≤ 45%

<sup>§</sup> For additional particle size information, please refer to the [Particle Size Distribution Cross Reference Chart](#) (Form No. 45-D00954-en).

<sup>1</sup> Contractual value


<sup>2</sup> Appears in current USP/NF.

DuPont™ Duolite™ AP143/1073 complies with the compendial specifications for Cholestyramine Resin USP when tested in conformance to the compendial test methods presented in current USP/NF. The resin is described as a “White to buff colored, hygroscopic, fine powder. Is odorless or has not more than a slight amine-like odor. Insoluble in water, in alcohol, in chloroform and in ether.”

## Identification

DuPont™ Duolite™ AP143/1073 can be identified by infrared spectroscopy, as shown in the example Figure 1.


Figure 1: Duolite™ AP143/1073 IR Spectrum

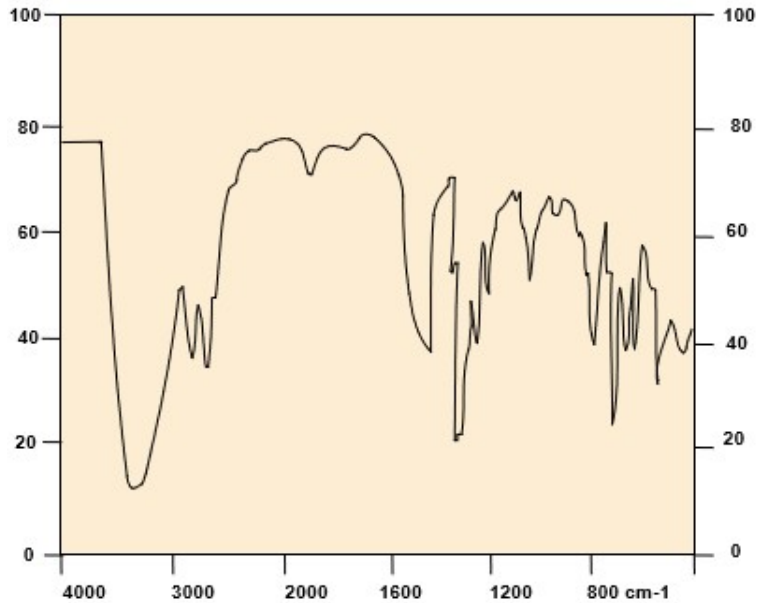

## Chemical Properties

Duolite™ AP143/1073 is derived from a copolymer of styrene and divinylbenzene with quaternary ammonium functionality. The mobile, or exchangeable, anion is chloride which can be exchanged for, or replaced by, virtually any anionic species. Duolite™ AP143/1073 Resin is an insoluble salt of a strong base and a strong acid; hence, its ability to exchange anions is largely independent of pH. The chemical structure of Duolite™ AP143/1073 is shown below in Figure 2.

Figure 2: Chemical Structure

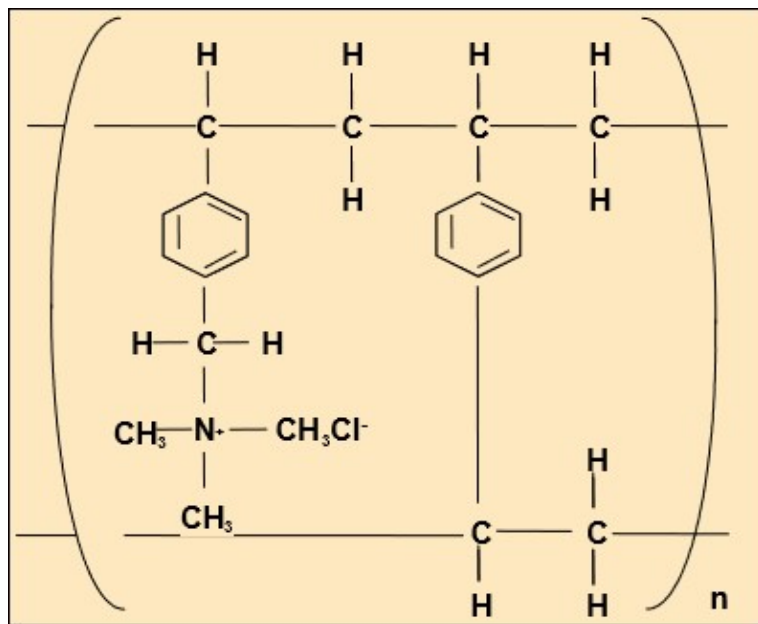

## **Applications**

- Taste Masking
- Drug Stabilization
- Controlled Release
- Active Ingredient

When used as a drug carrier, Duolite™ AP143/1073 provides a means for binding medicinal agents onto an insoluble polymeric matrix; this can be an effective technique to minimize taste and odor problems associated with the drug.

Controlled or sustained release properties can also be imparted to oral dosage formulations through the formation of resin-drug complexes (drug resins). The drug is released from the resin in vivo as the drug resin reaches equilibrium with the high electrolyte concentrations typically found in the gastrointestinal tract.

When used as an active ingredient, Duolite™ AP143/1073 resin binds bile acids; this leads to replenishment of the bile acids through increased catabolism of serum cholesterol, resulting in lowered serum cholesterol levels.

## **Cholestyramine Applications Reference List**

Irwin, W. J, R. MacHale, and P. J. Watts. (1990) Drug-delivery by ion exchange. Part VII: Release of acidic drugs from anionic exchange resin complexes. *Drug. Dev. Ind. Pharm.* 16(6):883-898

Sriwongjanya, Mongkol; Bodmeier, Roland; Effect of ion exchange resins on the drug release from matrix tablets; College Pharmacy, Freie Universitaet Berlin, Berlin, D-12169, Germany; *Eur. J. Pharm. Biopharm.* (1998), 46(3), 321-327

Polli, Gerald P. and Shoop, Clyde E., (Merck and Co. USA), 1976. Palatable cholestyramine coacervate compositions. Patent US 3,974,272.

Brauns H. A., Polli, Gerald P and Shoop, Clyde E., (Merck and Co., USA), 1974. Cholestyramine containing coacervate. *Ger. Offen DE 2,344,090.*

Kunin, Robert; Blood cholesterol reducing pharmaceutical composition; 1998; Patent US 5840339.

## **Product Stewardship**

DuPont has a fundamental concern for all who make, distribute, and use its products, and for the environment in which we live. This concern is the basis for our product stewardship philosophy by which we assess the safety, health, and environmental information on our products and then take appropriate steps to protect employee and public health and our environment. The success of our product stewardship program rests with each and every individual involved with DuPont products—from the initial concept and research, to manufacture, use, sale, disposal, and recycle of each product.

## Customer Notice

DuPont strongly encourages its customers to review both their manufacturing processes and their applications of DuPont products from the standpoint of human health and environmental quality to ensure that DuPont products are not used in ways for which they are not intended or tested. DuPont personnel are available to answer your questions and to provide reasonable technical support. DuPont product literature, including safety data sheets, should be consulted prior to use of DuPont products. Current safety data sheets are available from DuPont.

Please be aware of the following:

- **WARNING:** Oxidizing agents such as nitric acid attack organic ion exchange resins under certain conditions. This could lead to anything from slight resin degradation to a violent exothermic reaction (explosion). Before using strong oxidizing agents, consult sources knowledgeable in handling such materials.

Have a question? Contact us at:

[www.dupont.com/water/contact-us](http://www.dupont.com/water/contact-us)

All information set forth herein is for informational purposes only. This information is general information and may differ from that based on actual conditions. Customer is responsible for determining whether products and the information in this document are appropriate for Customer's use and for ensuring that Customer's workplace and disposal practices are in compliance with applicable laws and other government enactments. The product shown in this literature may not be available for sale and/or available in all geographies where DuPont is represented. The claims made may not have been approved for use in all countries. Please note that physical properties may vary depending on certain conditions and while operating conditions stated in this document are intended to lengthen product lifespan and/or improve product performance, it will ultimately depend on actual circumstances and is in no event a guarantee of achieving any specific results. DuPont assumes no obligation or liability for the information in this document. References to "DuPont" or the "Company" mean the DuPont legal entity selling the products to Customer unless otherwise expressly noted. NO WARRANTIES ARE GIVEN; ALL IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE ARE EXPRESSLY EXCLUDED. No freedom from infringement of any patent or trademark owned by DuPont or others is to be inferred.

© 2023 DuPont. DuPont™, the DuPont Oval Logo, and all trademarks and service marks denoted with ™, SM or ® are owned by affiliates of DuPont de Nemours Inc., unless otherwise noted.

